

BATERIE 675 HPX

DATE TEHNICE

SPECIFICATIE TEHNICA:

Marca/Produs de: Rayovac S.U.A.

Produs in: Marea Britanie

Tip: Zn-Aer , baterii speciale pentru aparate auditive, asigurand un sunet aproape fara distorsiuni

Gama: Extra (fosta Extra Advanced)

Tehnologii: Active Core Plus si Sound Fusion Technology

Cod produs: 675ZA, 675, PR44, B6754, DA675N6, s675, p675, B900PA, AC675E, HA675, 675AU, 7003ZD

Tensiune: 1.45 V

Capacitate: 700 mAh

Dimensiuni: maxim 5.4 x 11.6 mm (h x Ø)

Impedanta/rezistenta interna: <4 Ohm @ AC 1000Hz

Greutate: 1.79 grame/baterie

Durata standard de utilizare: 107 ore – 5mA background, 15mA 100ms pulsuri la fiecare 2h, 12 h/zi pana la 1.05V la 21°C si 50% R.H. (umiditate)


Durata conform IEC/ANSI High Drain: 57 ore – 8mA background, 24mA 100ms pulsuri o data la fiecare 2h, 12 h/zi pana la 1.05V la 21°C si 50% R.H. (umiditate)

Durata standard de depozitare: 4 ani in stare sigilata, la 21°C (dupa desigilare incep sa se auto-descarce)

Activare si utilizare: Inlaturati folia de plastic care acopera orificiile pentru aer. Inainte de a introduce bateria in aparat lasati-o un minut in aer liber pentru a-i asigura o activare corespunzatoare. Compartimentul bateriei trebuie sa permita o oxigenare adecvata



DIMENSIONS SHOWN ARE IEC
STANDARDS FOR ENGINEERING
TOLERANCES CONTACT
RAYOVAC.

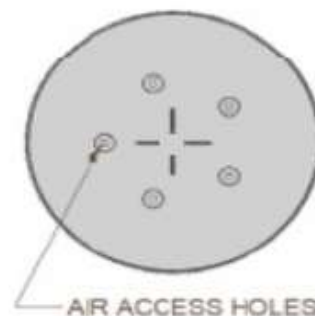

Standarde aplicabile: ANSI/NEDA - 7003ZD, IEC - PR44 / IEC 60086-1, IEC 60086-2 / ANSI C18.1M Parts 1 & 2

Ambalare: 6 buc/blister, 10 blistere/cutie (pretul este pentru o cutie cu 60 de baterii)

Protectia mediului: Bateria nu contine Mercur! Bateriile uzate nu se arunca la cosul de gunoi, ci se depun in urnele de colectare baterii uzate din marile magazine.

Typical IEC/ANSI Standard Discharge Performance

